




ELLB

Distribution Board for Resonant Piezo Motor Products

Operating Manual

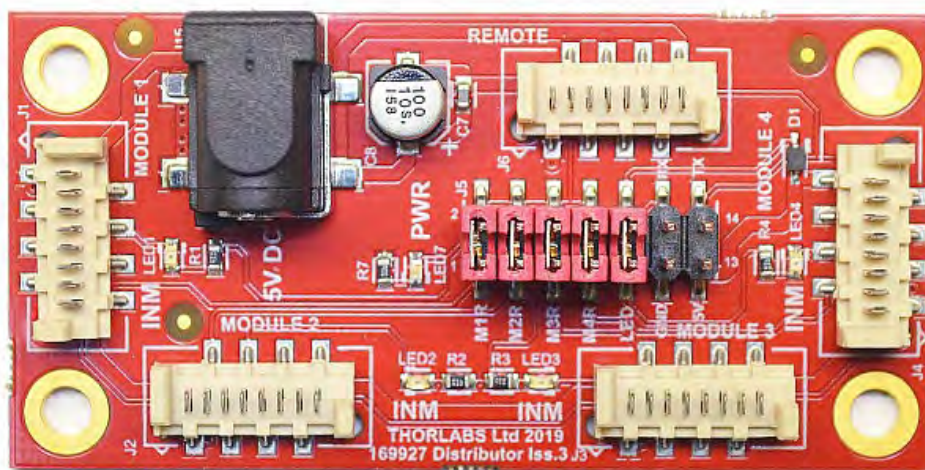

Original Instructions

Table of Contents

Chapter 1	Introduction	1
Chapter 2	Safety.....	2
2.1.	<i>General Warnings and Cautions</i>	2
Chapter 3	Operation	3
3.1.	<i>Environmental Conditions</i>	3
3.2.	<i>Getting Started</i>	3
3.2.1.	Distribution Board Description	4
3.2.2.	Interface Board Description.....	4
3.2.3.	Power Supply Options.....	5
3.2.4.	Connecting Devices	5
3.2.5.	Stage Connector Pin Out	6
3.3.	<i>Controlling the Stage</i>	7
3.3.1.	Software Control	7
3.3.2.	Communications Protocol	8
Chapter 4	Specifications.....	9
Chapter 5	Regulatory	10
5.1.	<i>Declarations of Conformity</i>	10
5.1.1.	For Customers in Europe	10
5.1.2.	For Customers In The USA.....	10
Chapter 6	Thorlabs Worldwide Contacts.....	11

Chapter 1 Introduction

The ELLB Distribution Board is part of the Thorlabs series of resonant piezo motor circuits and bare modules for OEM applications. The resonant piezo design of these motors offers fast response times and precise positioning, and are therefore particularly useful in scanning applications.

The high-speed digital signal processing (DSP) architecture supports a multi-drop serial communication protocol, and a set of digital IO lines allows the user to control the movement and state manually by switching the lines high (5V) or low (0V).

The board is designed to allow up to four elliptec resonant piezo motor devices to be controlled remotely via PC-based software, from a single GUI panel. The software can be downloaded from www.thorlabs.com. A compatible USB driver and source code are included in the software download package.

Connection to the host PC is via the USB interface board included with the Elliptec Evaluation Kits and power is supplied via an external 5V power supply.

Note

The module must be powered via an external 5V power supply unit (PSU). A suitable PSU is supplied in the ELL evaluation kits; the TPS101 T-Cube Power Supply is also compatible. For customers wanting to use a third party PSU, the Connector on the interface board is a DC Jack connector, 6.3mm OD (GND), 2.1mm ID (+5V).

Chapter 2 Safety

For the continuing safety of the operators of this equipment, and the protection of the equipment itself, the operator should take note of the Warnings, Cautions and Notes throughout this handbook and, where visible, on the product itself.

**Warning: Risk of Electrical Shock**

Given when there is a risk of electrical shock.

**Warning**

Given when there is a risk of injury to the user.

**Caution**

Given when there is a possibility of damage to the product.

Note

Clarification of an instruction or additional information.

2.1. General Warnings and Cautions

**Warning**

If this equipment is used in a manner not specified by the manufacturer, the protection provided by the equipment may be impaired. In particular, excessive moisture may impair operation.

The equipment is susceptible to damage from electrostatic discharge. When handling the device, anti-static precautions must be taken and suitable discharge appliances must be worn.



Spillage of fluid, such as sample solutions, should be avoided. If spillage does occur, clean up immediately using absorbent tissue.

Do not bend the PCB. A bending load in excess of 500 g applied to the board may cause the PCB to deform, which will degrade the performance of the controller.

During use do not place the PCB directly onto electro-conductive material e.g. an optical table top or breadboard.

Chapter 3 Operation

3.1. Environmental Conditions

**Warning**

Operation outside the following environmental limits may adversely affect operator safety.

Location	Indoor use only
Maximum altitude	2000 m
Temperature range	15°C to 40°C
Maximum Humidity	Less than 80% RH (non-condensing) at 31°C

To ensure reliable operation the unit should not be exposed to corrosive agents or excessive moisture, heat or dust.

If the unit has been stored at a low temperature or in an environment of high humidity, it must be allowed to reach ambient conditions before being powered up.

The unit is not designed to be used in explosive environments.

3.2. Getting Started

**Warning**

The safety of any system incorporating this equipment is the responsibility of the person performing the installation.

**Caution**

Although the module can tolerate up to 8 kV of air discharge, it must be treated as ESD sensitive device. When handling the device, anti-static precautions must be taken and suitable discharge appliances must be worn.



When power is applied, do not connect or disconnect the ribbon cable connecting the interface board to the distribution board. Always remove power before making connections.

The ribbon cable connector is made of plastic and is not particularly robust. Do not use force when making connections. Unnecessary or repeated plugging in and unplugging should be avoided or the connector may fail.

3.2.1. Distribution Board Description

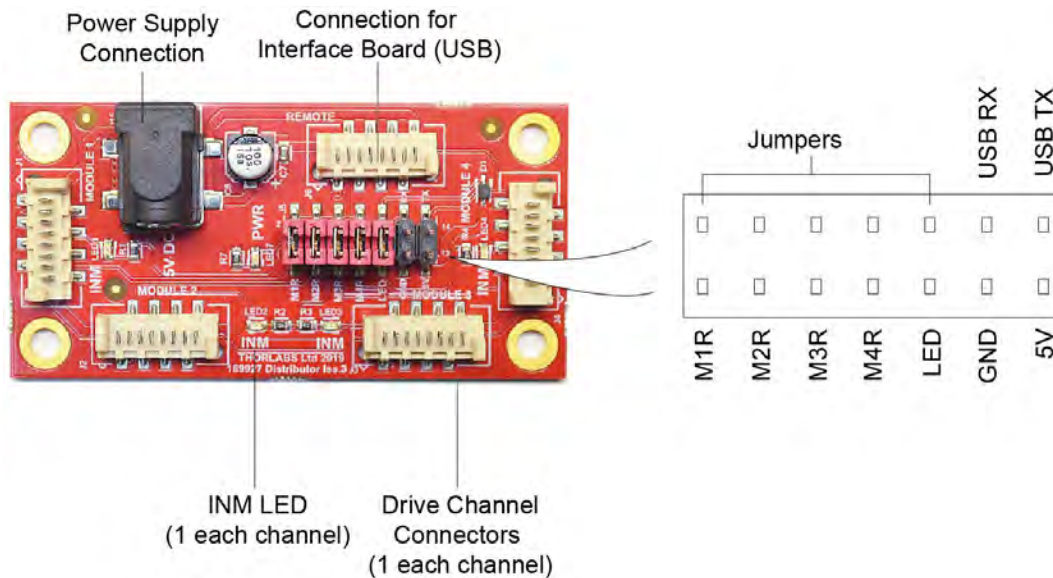


Figure 1 ELLB Distribution Board

Power Supply Connection – 5 V, 2 A, compatible with PSU supplied with ELLXXK kits or TPS101 T-Cube PSU. External diameter 6.3 mm, Internal diameter 2.1 mm.

Jumpers – M1R to M4R allow the reception of each of the four drive channels (M1R to M4R) to be enabled or disabled by the fitting or removal of the associated jumper. If the jumper is removed, the drive channel can still be addressed, but no response from the associated stage can be received.

LED allows the INM LEDs on the PCB to be enabled (jumper fitted) and disabled (jumper removed). If the jumper is fitted, the LED will be lit when the associated stage is moving.

USB RX and USB TX – These pins can be connected to a Raspberry or Arduino board instead of the USB Interface board to control the ELL stage modules connected to the distribution board.

GND and 5V are for use with an external power supply, 5 V 2 A.

Drive Channel Connectors (1 to 4) – These connectors allow either an interface board or a stage module to be connected to the distribution board.

INM LED – Lit when the associated stage is moving or busy.

3.2.2. Interface Board Description

The interface board supplied with the Elliptec Stage Evaluation Kits provides for USB connection to the host PC.



Caution

Moves cannot be initiated using the FW, BW and JOG buttons on the interface board. Moves must be performed using the Elliptec software running on a PC, or by writing a custom application using the messages described in the communications protocol document.

The PWR LED (LED1) is lit green when power is applied to the unit. The INM LED (LED2) is lit red when the device being driven is in motion.

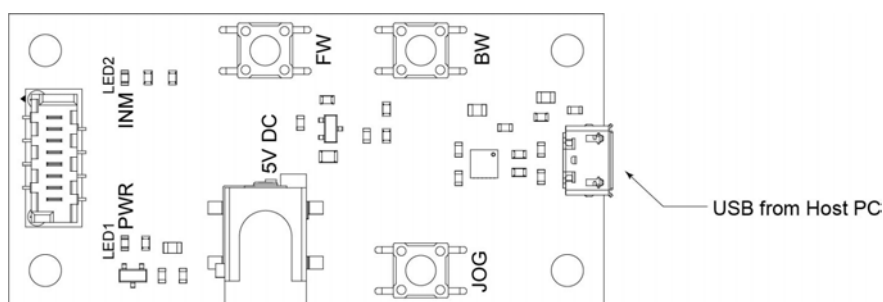


Figure 2 Handset Details

3.2.3. Power Supply Options

The distribution board is compatible with the TPS101 PSU available from www.thorlabs.com and with the power supply shipped with the ELLxx Evaluation Kits. These PSUs should be connected directly to the distribution board. This allows simultaneous motion of up to 2 modules.

The distribution board also features GND and 5V terminals (see section 3.2.1) for use with an external 5 V 2 A power supply like those available in most labs. If more current is needed (e.g. to move more than 2 devices at the same time) then a 5V 4A external power supply is available from Farnell (Part No 2815881).

**Caution**

As more devices are connected, simultaneous control of the units will require more current to be provided by the power supply.

Although the ELL6 can be operated standalone via USB power, when it is used with the ELLb then an external power supply is required.

3.2.4. Connecting Devices

When a device is first connected to the PC, it is assigned the default address '0'. The software can run multiple devices; however, before more than one device can be recognised, each device must be assigned a unique address. See below for a brief overview; detailed instructions are contained in the help file supplied with the software.

1. Switch ON the PC.
2. Connect an interface board (supplied with the individual ELL stage kits) to the distribution board.

**Caution**


The unit is easily damaged by connections with incorrect polarity. Pin 1 of the connector on the PCB is marked with an arrow (see Figure 4 and Figure 5.) which should be adjacent to the red wire in the connecting cable.

3. Connect the interface board to your PC.
4. Connect the first device to the distribution board.
5. Connect the distribution board to a 5V supply and switch 'ON'. (A 5V PSU is supplied with the ELL series stage kits, the TPS101 T-Cube supply is also compatible). Wait for the drivers to be installed
6. Run the Elliptec software and load the device.
7. Change the address of the first device.
8. Save the address of the first device.
9. Connect the next device to the distribution board.
10. In the ELL software GUI, press the "Disconnect" and "Connect" buttons to allow the software to see the new module.
11. Repeat steps 6 to 9 for the new module just connected.

Multiple devices can be controlled individually, either via the Elliptec software or by a third party application written using the messages detailed in the protocol document. Control via the handset is not possible when using the distribution board.

When powering down, remove power from the distribution board, before disconnecting the USB from the PC.

3.2.5. Stage Connector Pin Out

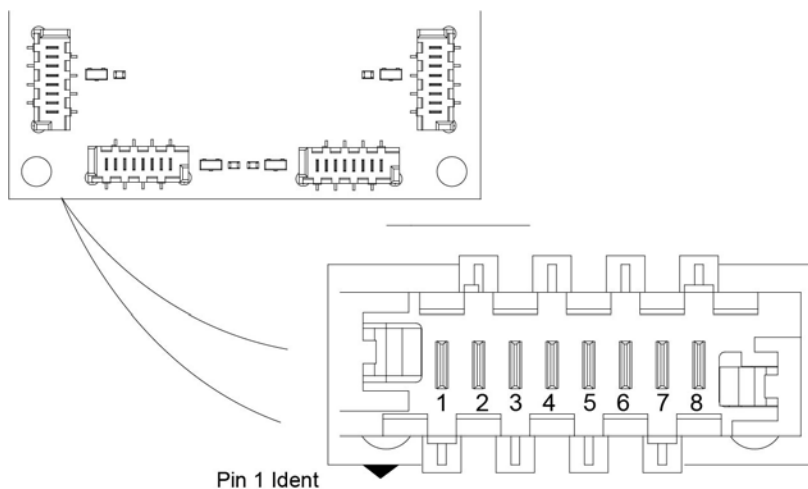

PIN	TYPE	FUNCTION
1	PWR	Ground
2	OUT	ODTX open drain transmit 3.3V TTL RS232
3	IN	RX receive 3.3V TTL RS232
4	OUT	In Motion, open drain active low max 5mA
5	IN	Not Connected
6	IN	Not Connected
7	IN	Not Connected
8	PWR	VCC +5V +/-10% 800mA (per module connected)

Connector model number MOLEX 90814-0808 Farnell order code 1518211

Mating connector model number MOLEX 90327-0308 Farnell order code 673160

Figure 3 Stage Connector Pin Out Details

3.3. Controlling the Stage

The stage can be controlled by the Elliptec software running on a PC, or by writing a custom application using the messages described in the communications protocol document. The modes of control are described in the following sections.

**Caution**

Moves cannot be initiated using the FW, BW and JOG buttons on the interface board.


3.3.1. Software Control

**Caution**


If multiple modules are connected, ensure that the address of each module has been changed (see section 3.2.4) before using the software.

When connected to the host PC, the stage can be controlled remotely, via the Elliptec software.

1. Download the Elliptec software from the downloads section at www.thorlabs.com. Double click the saved .exe file and follow the on-screen instructions.
2. Connect the hand-held controller to the stage unit.
3. Connect the hand-held controller to the 5V Power Supply and switch on.
4. Connect the hand-held controller to the PC USB port and wait for the drivers to be installed.
5. Run the Elliptec software.

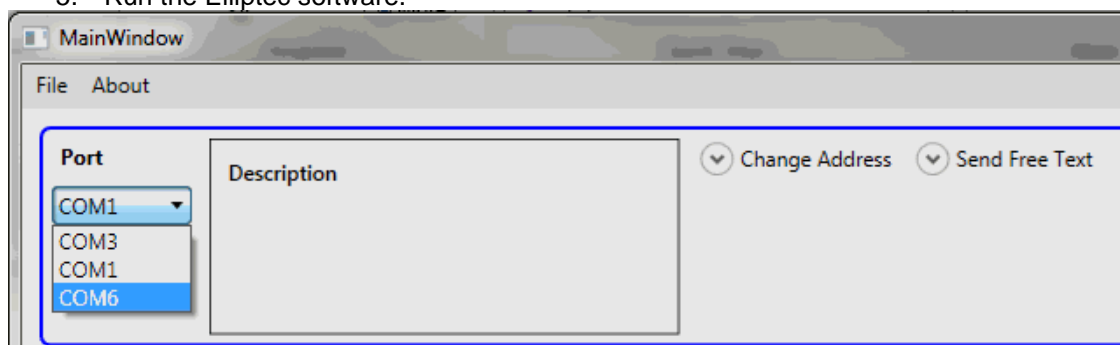

6. In the top left of the GUI panel displayed, select the COM port to which the device is connected, and click 'Connect'. The software will search the comms bus and enumerate the device.
7. Click the Home Offset 'Get' button, then click 'Home' to home the stage.
8. The GUI and device are now ready for use. See the helpfile supplied with the software for more information.

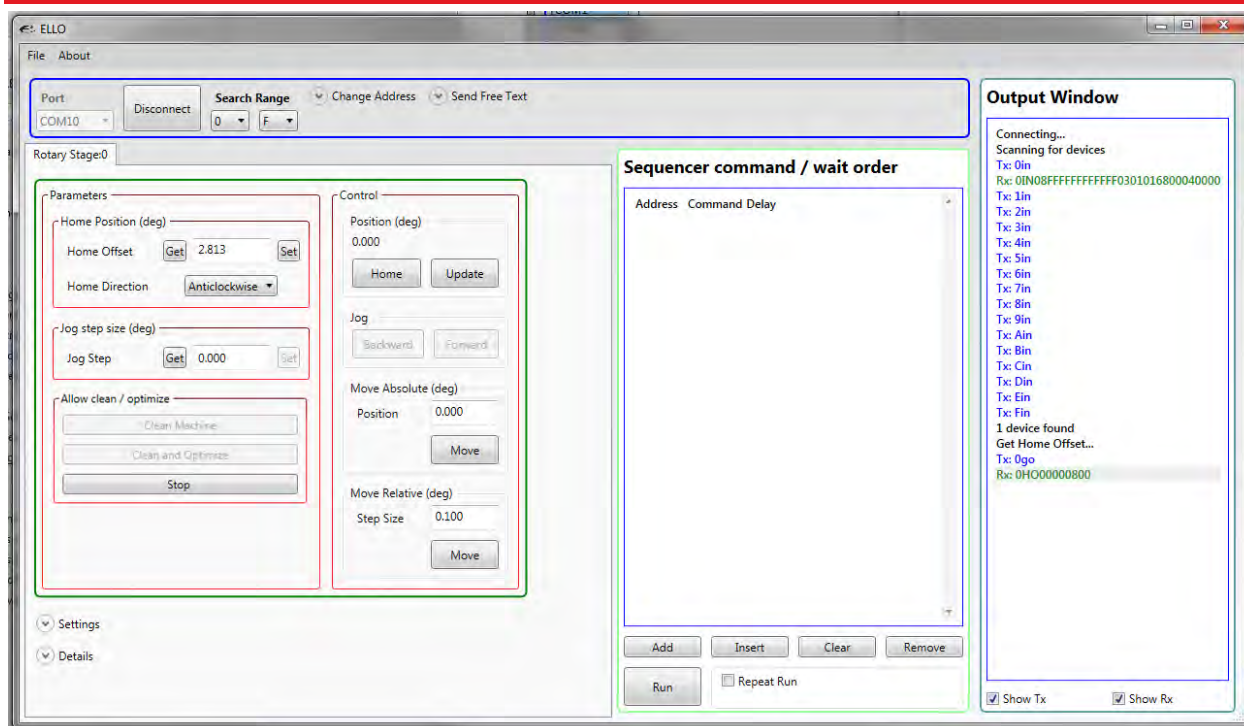

3.3.2. Communications Protocol

Custom move applications can be written in languages such as C# and C++.

The communication bus allows multi-drop communication with speeds at 9600 baud, 8 bit data length, 1 stop bit, no parity.

Protocol data is sent in ASCII HEX format, while module addresses and commands are mnemonic character (no package length is sent). Modules are addressable (default address is "0") and addresses can be changed and/or saved using a set of commands. Lower case commands are sent by user while upper case commands are replies by the module.

Please refer to the communications protocol manual for more detail about commands and data packet formats.


Chapter 4 Specifications

General Specifications	
Rated Voltage	4.5 to 5.5 V
Typical Current Consumption During Movement (per module connected)	800 mA
Max Current (of board)	4.0 A
Standby Current (per module connected)	0.05 A
8-Conductor Ribbon Cable Length (Supplied)	250 mm
8-Conductor Ribbon Cable Length (Maximum)	500 mm
Operating Temperature Range	15 to 40 °C (59 to 104 °F)
Weight	0.011 kg

Chapter 5 Regulatory

5.1. Declarations of Conformity

5.1.1. For Customers in Europe

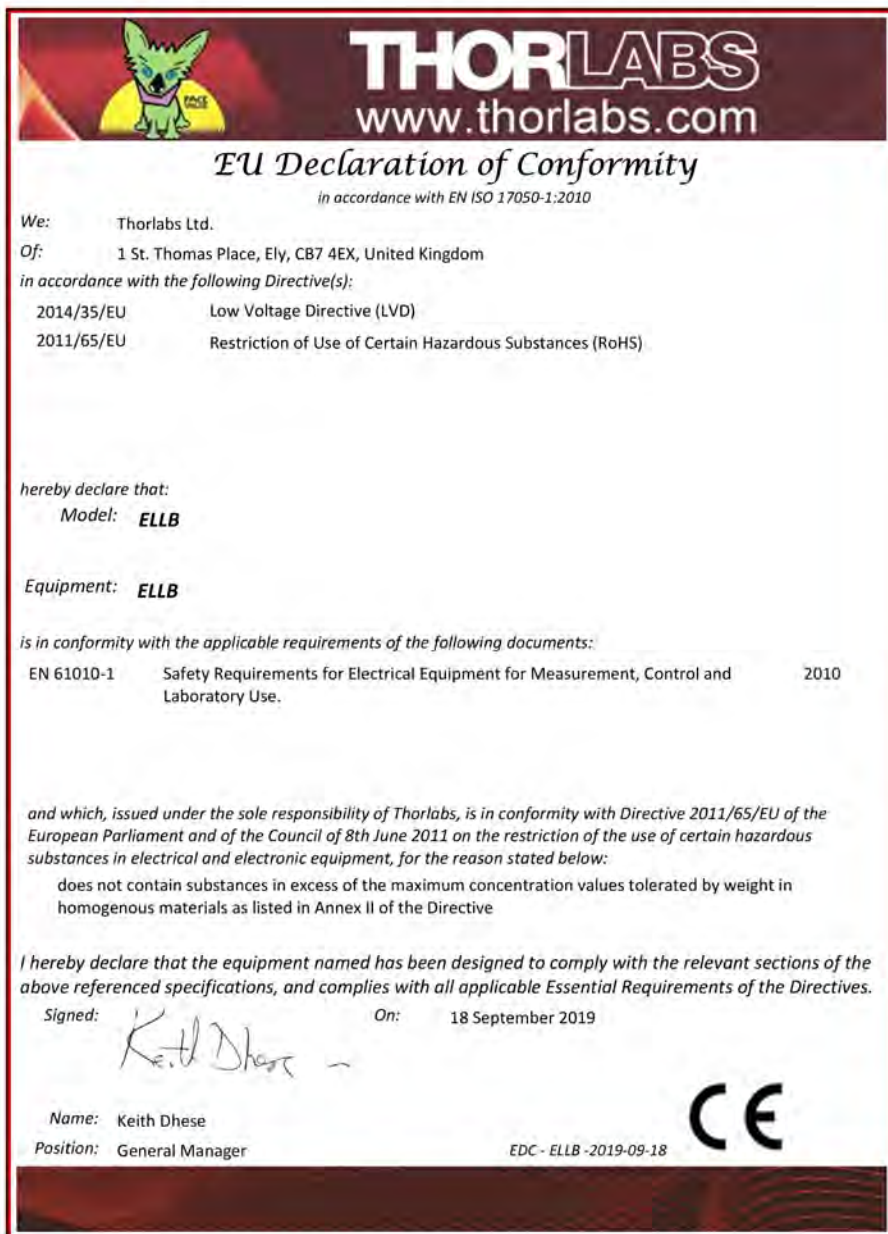

THORLABS
www.thorlabs.com

EU Declaration of Conformity
in accordance with EN ISO 17050-1:2010

We: Thorlabs Ltd.
Of: 1 St. Thomas Place, Ely, CB7 4EX, United Kingdom
in accordance with the following Directive(s):

2014/35/EU	Low Voltage Directive (LVD)
2011/65/EU	Restriction of Use of Certain Hazardous Substances (RoHS)

hereby declare that:
Model: **ELLB**

Equipment: **ELLB**


is in conformity with the applicable requirements of the following documents:

EN 61010-1	Safety Requirements for Electrical Equipment for Measurement, Control and Laboratory Use.	2010
------------	-------------------------------------------------------------------------------------------	------

and which, issued under the sole responsibility of Thorlabs, is in conformity with Directive 2011/65/EU of the European Parliament and of the Council of 8th June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, for the reason stated below:

does not contain substances in excess of the maximum concentration values tolerated by weight in homogenous materials as listed in Annex II of the Directive

I hereby declare that the equipment named has been designed to comply with the relevant sections of the above referenced specifications, and complies with all applicable Essential Requirements of the Directives.

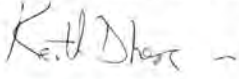Signed:  On: 18 September 2019

Name: Keith Dhese
Position: General Manager

EDC - ELLB -2019-09-18

CE

5.1.2. For Customers In The USA

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

Changes or modifications not expressly approved by the company could void the user's authority to operate the equipment.

Chapter 6 Thorlabs Worldwide Contacts

For technical support or sales inquiries, please visit us at www.thorlabs.com/contact for our most up-to-date contact information.



USA, Canada, and South America

Thorlabs, Inc.
sales@thorlabs.com
techsupport@thorlabs.com

Europe

Thorlabs GmbH
europe@thorlabs.com

France

Thorlabs SAS
sales.fr@thorlabs.com

Japan

Thorlabs Japan, Inc.
sales@thorlabs.jp

UK and Ireland

Thorlabs Ltd.
sales.uk@thorlabs.com
techsupport.uk@thorlabs.com

Scandinavia

Thorlabs Sweden AB
scandinavia@thorlabs.com


Brazil

Thorlabs Vendas de Fotônicos Ltda.
brasil@thorlabs.com

China

Thorlabs China
chinasales@thorlabs.com

Thorlabs verifies our compliance with the WEEE (Waste Electrical and Electronic Equipment) directive of the European Community and the corresponding national laws. Accordingly, all end users in the EC may return “end of life” Annex I category electrical and electronic equipment sold after August 13, 2005 to Thorlabs, without incurring disposal charges. Eligible units are marked with the crossed out “wheelie bin” logo (see right), were sold to and are currently owned by a company or institute within the EC, and are not disassembled or contaminated. Contact Thorlabs for more information. Waste treatment is your own responsibility. “End of life” units must be returned to Thorlabs or handed to a company specializing in waste recovery. Do not dispose of the unit in a litter bin or at a public waste disposal site.

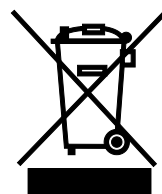

Annex I



THORLABS
www.thorlabs.com

